


Product Code . EL-BLE-10545

Ganong Potometer Apparatus

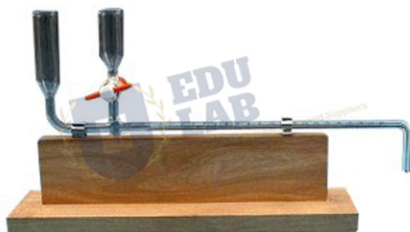

Description

Ganong Potometer Apparatus

For demonstration of transpiration rates and rates of water absorption by a cut shoot under various external conditions.

The open end of the capillary tube is dipped into a beaker of water.




A plant shoot is sealed into a stopper which fits the end tube and the whole apparatus is filled with water.

The position of the air bubble may easily be adjusted by opening the glass stopcock to admit water from the reservoir.

The whole apparatus is mounted on a varnished wooden base by spring clips.

Base: 20 x 80mm.

Capillary Bore: 1mm.

Educational Lab Equipments,
#449, HSIIDC, Industrial Area, Saha, Haryana
Direct Contact Details  +91-98173-19615  sales@educational-equipments.com
 www.educational-equipments.com